

AccuDiag™ Chikungunya IgG ELISA Kit

REF 8112-P1

IVD See External Label 2°C 96 Tests

Chikungunya IgG ELISA	
Principle	Indirect ELISA
Detection	Qualitative
Sample	4 µL serum/plasma
Incubation Time	100 minutes
Sensitivity	100%
Specificity	100%
Shelf Life	12 Months from the manufacturing date

PRODUCT FEATURES

- Very easy to use with little training
- Highly specific and consistent Assay
- Provides accurate results quickly
- Reading of results both visually and as absorbance data

INTENDED USE

The Diagnostic Automation, Inc. (DAI) Chikungunya IgG ELISA is designed for the qualitative detection of IgG antibodies present in human serum targeting chikungunya virus E2/E1 proteins. This test is for the presumptive clinical laboratory diagnosis of chikungunya virus infection. This assay is intended for use only in patients with clinical symptoms consistent with chikungunya infection. Positive results should be confirmed by following the guidelines recommended by the CDC for this disease.

SIGNIFICANCE AND SUMMARY

Chikungunya virus (CHIKV) is a mosquito-borne alphavirus first isolated during an epidemic in Tanzania in 1953 [1] [2]. Acute symptoms of the disease are

easily confused with dengue fever, and include sudden onset of high fever and polyarthralgia. The arthralgia can be debilitating and may last for periods of months or even years after the initial infection [3] [4]. As chikungunya may co-circulate with other diseases spread by *Aedes spp.* mosquitoes (e.g., dengue and yellow fever), it may be difficult to quickly and accurately diagnose the disease. CHIKV has re-emerged in recent years in a number of countries, causing outbreaks in Kenya, India, Italy and the Americas, resulting in significant healthcare consequences. Outbreaks on island populations (Lamu, Union of the Comoros) have resulted in epidemics that have infected greater than 50% of the population [3].

The DAI Chikungunya IgG ELISA is an enzyme linked sandwich-type immunoassay for the detection of human IgG antibodies targeting the CHIKV E2/E1 envelope glycoproteins. Polystyrene microtiter wells are pre-coated with capture antibodies for human IgG. Positive Control, Cut-Off Control, Negative Control and unknown test samples are diluted into a sample dilution buffer and then added to the ELISA plate. After incubation and washing, a subsequent ready-to-use CHIKV Antigen [5] [6] is added to each well. After a subsequent incubation and wash step, an HRP conjugated monoclonal antibody specific for CHIKV E2/E1 complex [7] is added to each well. After washing, wells are incubated with a tetramethylbenzidine (TMB) substrate. An acidic Stop Solution is then added and the degree of enzymatic turnover is determined by the absorbance (optical density) measurement at 450 nanometers. If human IgG antibodies targeting the CHIKV envelope glycoproteins are present, a complex is formed consisting of the IgG, antigen, and conjugate. If IgG antibodies targeting the CHIKV envelope glycoproteins are not present, then the complex of the antigen and conjugate are washed away.

SPECIMEN COLLECTION & PREPARATION

- Human serum must be used with this assay. Reagents have not been optimized, or tested with whole blood or plasma so these blood forms cannot be tested directly.
- Remove serum from the clot of red blood cells as soon as possible to avoid hemolysis.
- Testing should be performed as soon as possible after collection. Do not leave sera at room temperature for prolonged periods.
- Serum should be used and the usual precautions for venipuncture should be observed. The samples may be stored at 2-8°C for up to 7 days or frozen at -20°C or lower for up to 30 days. To maintain long-term longevity of the serum, store at -70°C. Avoid repeated freezing and thawing of samples.
- Frozen samples should be thawed to room temperature and mixed thoroughly by gentle swirling or inversion prior to use. Always quick-spin before use.
- If sera are to be shipped, they should be packed in compliance with Federal Regulations covering transportation of infectious agents.
- Do not use sera if any indication of growth is observed.

MATERIALS AND COMPONENTS

Materials provided with the kit
The DAI Chikungunya IgG ELISA contains sufficient reagents for one plate of 96 wells (12 x 8 strips) for human IgG targeting CHIKV. This is sufficient for testing a maximum of 90 unknown samples for human IgG, with controls included in duplicate.



Warning: Do not use any reagents where damage to the packaging has occurred.

- Coated Microtiter Test Strips for IgG (1 plate containing 12, 1x8 strips for human IgG):** ELISA plate strip holder with 96 (12x8 strips) coated polystyrene microtiter wells. Store at 2-8°C until expiry.
- Chikungunya IgG Negative Control (1x50µL):** The negative control will aid in verifying the validity of the kit. Store at 2-8°C until expiry. Centrifuge briefly prior to use to sediment any precipitate.
- Chikungunya IgG Positive Control (1x50µL):** The positive control will aid in verifying the validity of the kit. Store at 2-8°C until expiry. Centrifuge briefly prior to use to sediment any precipitate.
- Chikungunya IgG Cut-Off Control (1x50µL):** The cut-off control will aid in verifying the validity of the kit and establishing the threshold for reactive samples. Store at 2-8°C until expiry. Centrifuge briefly prior to use to sediment any precipitate.
- Sample Dilution Buffer for Chikungunya (2x25mL):** This solution is used for diluting all serum samples and controls prior to testing in the ELISA. Store at 2-8°C until expiry.
- Chikungunya Antigen (1x9mL):** This vial contains ready-to-use chikungunya antigen that comprises the chikungunya envelope glycoproteins. This component must be removed from the kit box and stored frozen. Thaw to room temperature before use. After each use, the antigen must be stored at -20°C or below until expiry. Avoid repeated freeze-thaw cycles (no more than 5).
- Conjugate Diluent for Chikungunya (1x9mL):** This solution is used for adding the enzyme conjugate to the ELISA assay. Do not use without first adding the 100X Conjugate, as described in the **Preparation of Reagents** section. Store at 2-8°C until expiry.
- 100X Conjugate for Chikungunya (1x150µL):** This contains horseradish peroxidase-labeled monoclonal antibody targeting chikungunya virus. Mix well prior to use. The 100X Conjugate is added to the Conjugate Diluent before use. Store the undiluted 100X conjugate at 2-8°C until expiry.
- 10X Wash Buffer (1x120mL):** One bottle of 10x concentrated Wash Buffer to be used as directed in Assay Procedure. Store at 2-8°C until expiry.
- Liquid TMB Substrate (1x12mL):** To be used as directed in Assay Procedure. Store at 2-8°C until expiry.

Note: This substrate is light-sensitive and should be stored in the original bottle.

- Stop Solution (1x9mL):** To be used to terminate the reaction as directed in Assay Procedure. Store at 2-8°C until expiry.

Caution: This is a strong acid, wear protective gloves, mask and safety glasses. Dispose of all materials according to safety rules and regulations.

Materials required but not provided

- ELISA Spectrophotometer capable of absorbance measurement at 450 nm
- Biological or High-Grade Water
- Vacuum Pump
- Automatic Plate Washer
- 37°C Incubator without humidification or CO₂
- 1-10 µL Single-Channel Pipettors, 50-200 µL Single- and Multi-Channel Pipettors
- Polypropylene tubes or 96 well dilution plates
- Parafilm or similar plate cover
- Timer
- Vortex

REAGENT PREPARATION

Bring all kit reagents and specimens to room temperature (~25°C) before use. Thoroughly mix the reagents and samples before use by gentle inversion.

CAUTION: this kit has not been optimized by DAI for use with any particular automated ELISA processing system. Use with an automated Elisa processing system will require proper validation to ensure results are equivalent to the expectations described in this package insert. Modifications to the protocol of these systems and/or different volumes of reagents may be required.


Preparation of Reagents:

- Preparation of 1X Wash Buffer**
Dilute the 10X Wash Buffer to 1X using Biological or High-Grade Water. To prepare a 1X wash buffer solution, mix 120 mL 10X wash buffer with 1080 mL distilled (or deionized water). Mix thoroughly to ensure that any precipitate is dissolved and that the solution is uniform. Once diluted to 1X, the solution can be stored at room temperature for up to 6 months. Check for contamination prior to use. Discard if contamination is suspected.
- Microtitration Wells**
Select the corresponding strips depending upon the number of samples to test. Evaluate the plate to ensure the strips are flat and properly aligned to the plate. The remaining unused wells should be repackaged and sealed immediately with the supplied desiccant and stored at 2-8°C until ready to use or expiration.
- Preparation of the Chikungunya Ready To Use Enzyme Conjugate**
Add 90 µL of 100X Conjugate for Chikungunya directly to the 9 mL bottle of Conjugate Diluent for Chikungunya (1 part: 100 parts). Mix by inverting solution several times. For smaller volumes, aliquot the required volume of sample diluent into a separate, clean polypropylene test tube and add the appropriate volume of 100X Conjugate for Chikungunya to the diluent in the test tube. Prepare only the required amount needed for a given test run.

ASSAY PROCEDURE

- Positive, Cut-Off and Negative Controls should be assayed in duplicate (and run each time the assay is performed on every plate). Unknown serum samples may be assayed in singlet. However, it is recommended to run samples in duplicate until the operator is familiar with the assay.
- Dilute the Controls and unknown serum samples in the Sample Dilution Buffer for Chikungunya. Samples should be diluted 1/100 into a clean dilution plate well or separate clean polypropylene test tube. For example, dilute 4 µL of sample into 396 µL of Sample Dilution Buffer. Mix well by pipetting up and down several times.
- Add 50 µL of the 1/100 diluted Controls and unknown serum samples to the appropriate ELISA plate wells.
- Cover the plate with parafilm, as shown below.

Note: This is to ensure even temperature distribution in all wells from bottom and sides; any extra parafilm can be cut off once the top is sealed to block evaporation.

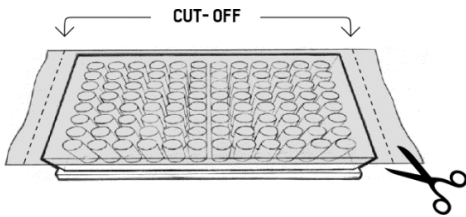


- Incubate at 37°C for 30 minutes (± 1 minutes) in an incubator. The incubator should be properly calibrated and verified to maintain 37°C ± 2°C using an external, reference thermometer.

Note: Do not stack plates on top of each other. They should be spread out as a single layer. This is very important for even temperature distribution. Do not use CO₂ or other gases. Do not place plates in contact with any wet substances such as wet paper towels etc.



CORRECT METHOD

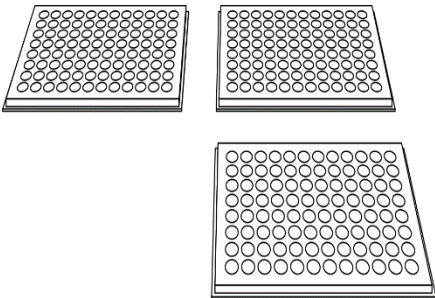

- After the incubation, wash the plate 6 times with an automatic plate washer using 1x Wash buffer. Use 300 µL per well in each wash cycle.
- Add 50 µL of the Chikungunya Antigen to each ELISA plate well.
- Cover the plate with parafilm or similar cover.
- Incubate at 37°C for 30 minutes (± 1 minute) in an incubator.
- During the incubation, prepare a fresh volume of Ready To Use Enzyme Conjugate as described in the Preparation of Reagents section (e.g., add 90 µL of 100X Conjugate for Chikungunya directly to the 9 mL bottle of Conjugate Diluent for Chikungunya, 1 part: 100 parts). Prepare only the minimum amount of Ready To Use Enzyme Conjugate necessary for the assay run.
- After the incubation, wash the plate 6 times with an automatic plate washer using 1x Wash buffer. Use 300 µL per well in each wash cycle.
- Add 50 µL of the Ready To Use Enzyme Conjugate to each ELISA plate well.
- Cover the plate with parafilm or similar cover.
- Incubate at 37°C for 30 minutes (± 1 minute) in an incubator.

- After the incubation, wash the plate 6 times with an automatic plate washer using 1x Wash buffer. Use 300 µL per well in each wash cycle.
- Apply 75 µL of Liquid TMB substrate into all wells using a multi-channel pipetter. **Note:** The Liquid TMB is light-sensitive. Avoid exposure to unnecessary light sources.
- Incubate the plate in the dark, at room temperature for 10 minutes (± 30 seconds).
- Add 50 µL of Stop Solution into all wells using a multi-channel pipetter and let the plate stand, uncovered at room temperature for 1 minute.
- Read the optical density value at 450nm (OD₄₅₀) with a Microplate reader. DO NOT SUBTRACT OR NORMALIZE ANY BLANK VALUES OR WELLS.
- Analyze the data as shown in the Quality Control section and Interpretation of Results section.

Procedure Summary

Dilute Controls and Samples 1/100 into the Sample Dilution Buffer (4µL:396µL). Controls should be assayed in duplicate.

Add 50µL of diluted Controls and Samples to the appropriate ELISA plate wells.

Cover with parafilm, incubate 30 minutes, 37°C.

Wash plate 6x. Add 50 µL of the Chikungunya Antigen to each ELISA plate well. Make sure to return the antigen to storage at -20°C or below.

Cover with parafilm, incubate 30 minutes, 37°C.

Prepare a fresh volume of Ready To Use Enzyme Conjugate (e.g., add 90 µL of 100X Conjugate for Chikungunya directly to the 9 mL bottle of Conjugate Diluent for Chikungunya). Prepare only the minimum volume necessary for the assay run.

Wash plate 6x. Add 50 µL of Ready To Use Enzyme Conjugate to each ELISA plate well.

Cover with parafilm, incubate 30 minutes, 37°C.

Wash plate 6x. Add 75 µL of TMB to all wells. Avoid exposure to light. Incubate 10 minutes at room temperature.

Add 50 µL of Stop solution to each well. Measure the OD₄₅₀ (do NOT subtract background or reference wavelengths).

RESULTS

Samples with ISR values ≥ 1.0 are considered “Reactive” and samples with ISR values < 1.0 are considered “Non-Reactive”. Any samples in the range of 0.90 < ISR Value < 1.10, should be considered retested and the sample evaluated in duplicate to verify the sample status.

ISR Value	Interpretation
< 1.0	“Non-Reactive”. Sample is considered to not have IgG antibodies specific for chikungunya
≥ 1.0	“Reactive”. Sample is considered to have IgG antibodies specific for chikungunya present
0.9 - 1.10	Repeat test sample in duplicate to verify results



QUALITY CONTROL

Each kit contains positive, cut-off and negative control sera. Positive, cut-off and negative controls must be run on each plate tested. Acceptable optical density (OD) values for these controls are found on the specification table below. The negative and positive controls are intended to monitor for substantial reagent failure. The test is invalid and must be repeated if either of the controls do not meet the specifications. The cut-off control is provided to establish a threshold above which a sample is considered reactive. Quality Control (QC) procedures must be performed in conformance with local, state, and/or federal regulations or accreditation requirements and the user's own laboratory's standard QC procedures. It is recommended that the user refer to CLSI C24-A and 42 CFR 493.1256 for guidance on appropriate QC practices. The results below are given strictly for guidance purposes only. Analysis is applicable when using RAW spectrophotometric readings only and where automatic subtraction of water or reagent blanks is not employed.

Calculation of the Chikungunya IgG Negative Control:

Calculate the mean (average) Negative Control value for the Chikungunya IgG test.

Example 1: Mean Negative Control

OD ₄₅₀	
Replicate 1	0.084
Replicate 2	0.094
Sum	0.178

Mean OD₄₅₀ = 0.178 ÷ 2 = 0.089

Calculation of the Chikungunya IgG Positive Control:

Calculate the mean (average) Positive Control value for the Chikungunya IgG test.

Example 2: Mean Positive Control

OD ₄₅₀	
Replicate 1	1.540
Replicate 2	1.378
Sum	2.918

Mean OD₄₅₀ = 2.918 ÷ 2 = 1.459

Calculation of the Chikungunya IgG Cut-off Control Value:

Calculate the mean (average) Cut-Off Control value for the Chikungunya IgG test.

Example 3: Mean Cut-Off Control

OD ₄₅₀	
Replicate 1	0.182
Replicate 2	0.194
Sum	0.376

Mean OD₄₅₀ = 0.376 ÷ 2 = 0.188
Cut-Off Control Value = 0.188

Calculation of the Chikungunya IgG Discrimination Capacity:

Calculate the Discrimination Capacity (R_{PC/NC}) for the assay. The Discrimination Capacity (R_{PC/NC}) for the ELISA is defined as the ratio of the Positive Control to the Negative Control.

Example 4: Discrimination Capacity (R_{PC/NC})

OD ₄₅₀	
Mean Positive OD ₄₅₀	1.459
Mean Negative OD ₄₅₀	0.089
Ratio	16.39

R_{PC/NC} = 1.459 ÷ 0.089 = 16.39

Calculation of the Chikungunya IgG Immune Status Ratio (ISR) Value:


Calculate the Immune Status Ratio (ISR) for each sample in the assay. The ISR value for each sample is defined as the ratio of the test sample OD₄₅₀ to the mean Cut-Off Control value.

Example 5: ISR Value

OD ₄₅₀	
Mean Test Sample OD ₄₅₀	0.876
Cut-Off Control Value	0.188

ISR = 0.876 ÷ 0.188 = 4.66

Quality Control Requirements: The following criteria must be fulfilled for the Diagnostic Automation, Inc. Chikungunya IgG ELISA in order for the assay run to be considered valid.

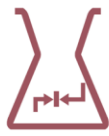

Quality Control Requirements	
	Requirement
Mean Negative OD ₄₅₀	≤ 0.200
Mean Positive OD ₄₅₀	≥ 0.500
Mean Cut-Off OD ₄₅₀	> Mean Neg. Control OD ₄₅₀
R _{PC/NC}	≥ 5.00

For accurate results do not subtract "blank" well values or reference wavelenghts

LIMITATIONS OF THE ASSAY

1. For in vitro diagnostic use only. Not for sale or distribution in the United States of America.
2. The possibility of false positive and false negative results must be considered.
3. Serological cross-reactivity with other alphavirus samples must be considered. Certain sera (e.g., those positive for Eastern equine encephalitis virus, human anti-mouse antibodies, etc.) may give false positive results. Therefore, any positive sera must be confirmed with other tests.

PRECAUTIONS

1. For professional *in vitro* diagnostic use only.
2. Not for sale or distribution in the United States of America.
3. All human source materials used in the preparation of controls have tested negative for antibodies to HIV 1&2, Hepatitis C and Hepatitis B surface antigen. However, no test method can ensure 100% efficiency or sensitivity. Therefore, all human controls and antigen should be handled as potentially infectious material. The Centers for Disease Control and Prevention and the National Institutes of Health recommend that potentially infectious agents be handled at the Biosafety Level 2.
4. A thorough understanding of this package insert is necessary for successful use of the product. Reliable results will only be obtained by using precise laboratory techniques and accurately following the package insert.
5. Do not mix various lots of any kit component within an individual assay.
6. Do not use any component beyond the expiration date shown on its label.
7. Avoid exposure of the reagents to excessive heat or direct sunlight during storage and incubation.
8. Some reagents may form a slight precipitate, mix gently before use.
9. Incomplete washing will adversely affect the outcome and assay precision.
10. To minimize potential assay drift due to variation in the substrate incubation time, care should be taken to add the stopping solution into the wells in the same order and speed used to add the TMB solution.
11. Avoid microbial contamination of reagents.

12. Avoid contamination of the TMB Substrate Solution with the Enzyme Conjugate-HRP.
13. Wear protective clothing, eye protection and disposable gloves while performing the assay. Wash hands thoroughly afterwards.
14. Use a clean disposable pipette tip for each reagent, Standard, Control or specimen.
15. Cover working area with disposable absorbent paper.

WARNING: Potentially biohazardous material

This kit may contain reagents made with human serum or plasma. The serum or plasma used has been heat-inactivated unless otherwise stated. Handle all sera and kits used as if they contain infectious agents. Observe established precautions against microbiological hazards while performing all procedures and follow the standard procedures for proper disposal of specimens.

CHEMICAL HAZARD:

Safety Data Sheets (SDSs) are available for all components of this kit. Review all appropriate SDS before performing this assay. Avoid all contact between hands and eyes or mucous membranes during testing. If contact does occur, consult the applicable SDS for appropriate treatment.


REFERENCES

1. S. C. Weaver, J. E. Osorio, J. A. Livengood, R. Chen and D. T. Stinchcomb, "Chikungunya virus and prospects for a vaccine," *Expert Rev. Vaccines*, vol. 11, no. 9, pp. 1087-1101, 2012.
2. B. L. Tang, "The cell biology of Chikungunya virus infection," *Cellular Microbiology*, vol. 14, no. 9, pp. 1354-1363, 2012.
3. J. E. Staples, R. F. Breiman and A. M. Powers, "Chikungunya Fever: An Epidemiological Review of a Re-Emerging Infectious Disease," *Emerging Infections*, vol. 49, pp. 942-948, 2009.
4. S. C. Weaver, "Arrival of Chikungunya Virus in the New World: Prospects for Spread and Impact on Public Health," *PLOST Neglected Tropical Diseases*, vol. 8, no. 6, p. e2921, 2014.
5. F. Nasar, G. Palacios, R. V. Gorchakov, H. Guzman, A. P. Travassos Da Rosa, N. Savji, V. L. Popov, M. B. Sherman, W. I. Lipkin, R. B. Tesh and S. C. Weaver, "Eilat virus, a unique alphavirus with host range restricted to insects by RNA replication," *PNAS*, vol. 109, no. 36, pp. 14622-14627, 2012
6. Nasar F, Gorchakov RV, Tesh RB, Weaver SC. Eilat Virus Host Range Restriction Is Present at Multiple Levels of the Virus Life Cycle. *J Virol*. 2015;89(2):1404-18.
7. Pal P, Dowd KA, Brien JD, Edeling MA, Gorlatov S, Johnson S, Lee I, Akahata W, Nabel GJ, Richter MK, Smit JM, Fremont DH, Pierson TC, Heise MT, Diamond MS. Development of a highly protective combination monoclonal antibody therapy against Chikungunya virus. *PLoS Pathog*. 2013;9(4):e1003312.




MANUFACTURER AND BRAND DETAILS

ISO 13485:2016



ISO 13485
Quality
Management for
Medical Devices
CERTIFIED

 Diagnostic Automation/Cortez Diagnostics, Inc.
21250 Califa Street, Suite 102 and 116,
Woodland Hills, California 91367 USA

Date Adopted	2023-09
Brand Name	AccuDiag™
REF 8112-P1	AccuDiag™ - Chikungunya IgG ELISA
EC REP	CEpartner4U, Esdoornlaan 13, 3951 DB Maarn, The Netherlands. www.cepartner4u.eu

Revision Date: 2019-09-24